


# Cotton Tipped Applicator

Wood handle

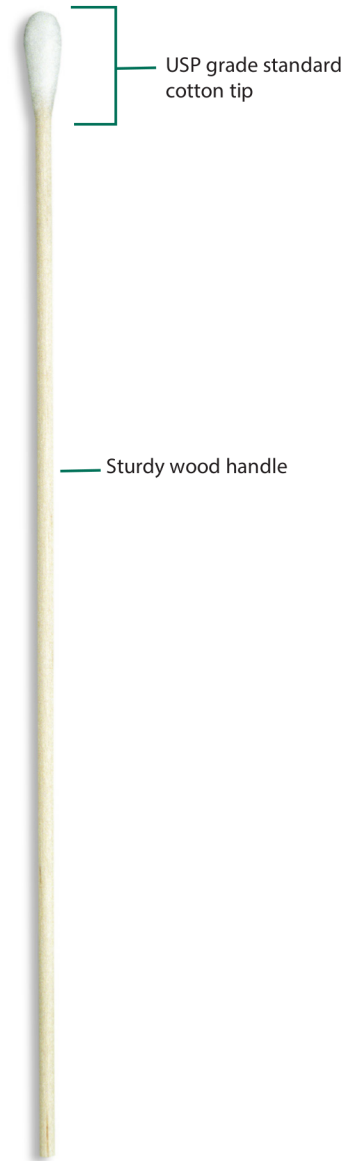

## Description

This absorbent cotton tipped applicator is produced with High Quality USP grade cotton which is carded and coiled onsite to maintain consistency. It is bonded securely to a rigid white birch handle by an aqueous based adhesive.

## Dimensions

6" Overall Length (76.2 mm)

### Tip

- 0.188" Dia. (4.78mm)
- 0.625" Length (15.88 mm)

### Handle

- 0.100" Dia. (2.54 mm)
- 5.812" Length (147.62 mm)

## Qualities

- Medical grade quality
- Absorbent tip
- Soft smooth tip
- Made in USA

## Applications

- Applying salves or ointments
- Wound care
- Application of topical dressings
- Specimen collection
- Applying adhesives or other liquids
- Removing dirt, contamination, fluxes or resins
- Cleaning with IPA or other solvents

Ref #	Packaging
806-WC Non-sterile	10/10/100 (10 boxes per case, 100 bags per box, 100 swabs per bag)
806-WCB Non-sterile	100/100 Bulk (100 bags per case, 100 swabs per bag)



**Puritan**

Puritan Medical Products Company LLC  
 Guilford, Maine 04443-0149 USA  
 ISO 9001 2008 ISO 13485:2003 CE

Tel 800 321-2313 (US & Canada Only)  
 or 207-876-3311

Web [puritanmedproducts.com](http://puritanmedproducts.com)  
 Email [sales@puritanmedproducts.com](mailto:sales@puritanmedproducts.com)

# Cotton Tipped Applicator

Wood handle

Ref #: 806-WC / 806-WCB

## NVR

NON-VOLATILE RESIDUE Test method: IDEMA method M7-98

NVR Solvent - Hexane

NVR 12 µg/swab

The minimum detection limit for NVR is 10 µg/swab

## ION CHROMATOGRAPHY TESTING RESULTS

IONIC CLEANLINESS

Test method: IDEMA M13-99 (anions) M12-99 (cations)

### Anions

Fluoride	ND
Chloride	0.9000 µg/swab
Nitrate	ND
Bromide	2.5000 µg/swab
Nitrite	0.3250 µg/swab
Sulfate	4.3000 µg/swab
Phosphate	11.2500 µg/swab


### Cations

Lithium	ND
Sodium	10.8000 µg/swab
Ammonium	0.7500 µg/swab
Potassium	13.7550 µg/swab
Magnesium	2.7450 µg/swab
Calcium	6.8250 µg/swab

ND - Non Detected

Test procedures available upon request.

## Distributed By:

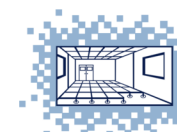

**HUTCHINS & HUTCHINS**  
 Your Cleanroom Supplier

Call for more info, pricing & availability | Toll Free: (800) 554-4736 | Local: (540) 949-6663