

**STX705P** Sterile Cotton Polystyrene Handle



**STX705W** Sterile Cotton Wood Handle



**STX763** Sterile Polyester



**STX764** Sterile Polyester Small



### Description

Texwipe's Sterile Cotton Swabs are constructed from 100% USP-grade cotton. Texwipe's Sterile Polyester Swabs are constructed from 100% USP-grade spun polyester. The swabs are bonded securely to a wooden or polystyrene shaft by an aqueous-based adhesive. All swabs are made to exacting and consistent tolerances using high-precision automated processes.

The easy-to-open, individually wrapped swabs ensure sterility at point-of-use. The swabs are visible through clear packaging which also offers visibility of seal integrity.

### Features & Benefits

- Gamma irradiated to a Sterility Assurance Level of 10<sup>-6</sup> in accordance with ANSI/AAMI/ISO 11137
- Lot code and expiration date clearly marked on each package and included in bar code for inventory control
- Certificates of Irradiation and Compliance attached to each case confirm irradiation dose met specifications
- Pre-sterilized swabs eliminate repackaging, documentation and validation of autoclaved swabs
- Packaging labeled for global use

Sterile swabs available upon request.

### Applications

- Diagnostic sampling
- Environmental monitoring
- Cleaning with solvents such as IPA

### Industries

- Pharmaceuticals
- Biologics
- Medical Device
- Microbiological Laboratories
- Diagnostic Testing Laboratories

### Products

See reverse for Physical Characteristics


### Physical Characteristics

	STX705P	STX705W	STX763	STX764
Head material	Cotton	Cotton	Polyester	Polyester
Head width	7.0 mm (0.276")	7.0 mm (0.276")	5.8 mm (0.228")	4.6 mm (0.181")
Head thickness	7.0 mm (0.276")	7.0 mm (0.276")	5.8 mm (0.228")	4.6 mm (0.181")
Head length	18.2 mm (0.717")	18.2 mm (0.717")	18.0 mm (0.709")	17.0 mm (0.669")
Handle material	Polystyrene	Wood	Polystyrene	Polystyrene
Handle width	2.5 mm (0.098")	2.5 mm (0.098")	2.5 mm (0.098")	2.5 mm (0.098")
Handle thickness	2.5 mm (0.098")	2.5 mm (0.098")	2.5 mm (0.098")	2.5 mm (0.098")
Handle length	134.0 mm (5.276")	134.0 mm (5.276")	134.0 mm (5.276")	134.0 mm (5.276")
Total swab length	155.0 mm (6.102")	155.0 mm (6.102")	155.0 mm (6.102")	155.0 mm (6.102")
Head bond	Adhesive	Adhesive	Adhesive	Adhesive
Handle color	White	Brown	White	White
Design notes	100% cotton head; long polystyrene handle	100% cotton head; long wooden handle	100% spun polyester head; long polystyrene handle	100% spun polyester head; long polystyrene handle

### Products

Number	Description	Packaging
<b>STX705P</b>	Sterile Cotton Swab with Polystyrene Handle	1 swab/sleeve; 50 swabs/box; 10 boxes/case
<b>STX705W</b>	Sterile Cotton Swab with Wood Handle	1 swab/sleeve; 50 swabs/box; 10 boxes/case
<b>STX705PT</b>	Sterile Dry Collection and Transport System, Cotton Swab with Polystyrene Handle	1 swab/tube; 100 tubes/box; 5 boxes/case
<b>STX763</b>	Sterile Polyester Swab	1 swab/sleeve; 50 swabs/box; 10 boxes/case
<b>STX764</b>	Sterile Small Polyester Swab	1 swab/sleeve; 50 swabs/box; 10 boxes/case
<b>STX764T</b>	Sterile Dry Collection and Transport System, Small Polyester Swab	1 swab/tube; 100 tubes/box; 5 boxes/case

### Distributed By:

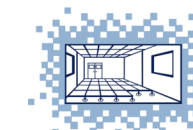

**HUTCHINS & HUTCHINS**<sup>™</sup>  
Your Cleanroom Supplier<sup>®</sup>

**Call for more info, pricing & availability | Toll Free: (800) 554-4736 | Local: (540) 949-6663**

**Texwipe Regulatory Certificates**  
Texwipe is registered: ISO 9001 • ISO 14001 • OHSAS 18001